

Display Format of the Selected Plots


When you add plots to the **Plots** area they are added next to each other in a **Gallery** plot view.

You can change the display to a family tree style hierarchy format by clicking the **Hierarchy** button in the **View** group of the **Analyze** tab. In this format the plots are displayed in a tree-like structure with a separate branch for each gate. Selecting the **Gallery** button will return you to the **Gallery** plot view.

A plot containing a region which is part of a gate equation is called a **Classifier** plot. A plot that has no gating regions on it is called an **Outcome** plot.

A pathway is drawn between each classifier plot and between other classifier plots and the subsequent gated outcome plots.

Ungated plots with no gating regions drawn on them appear outside the hierarchy as orphaned plots.

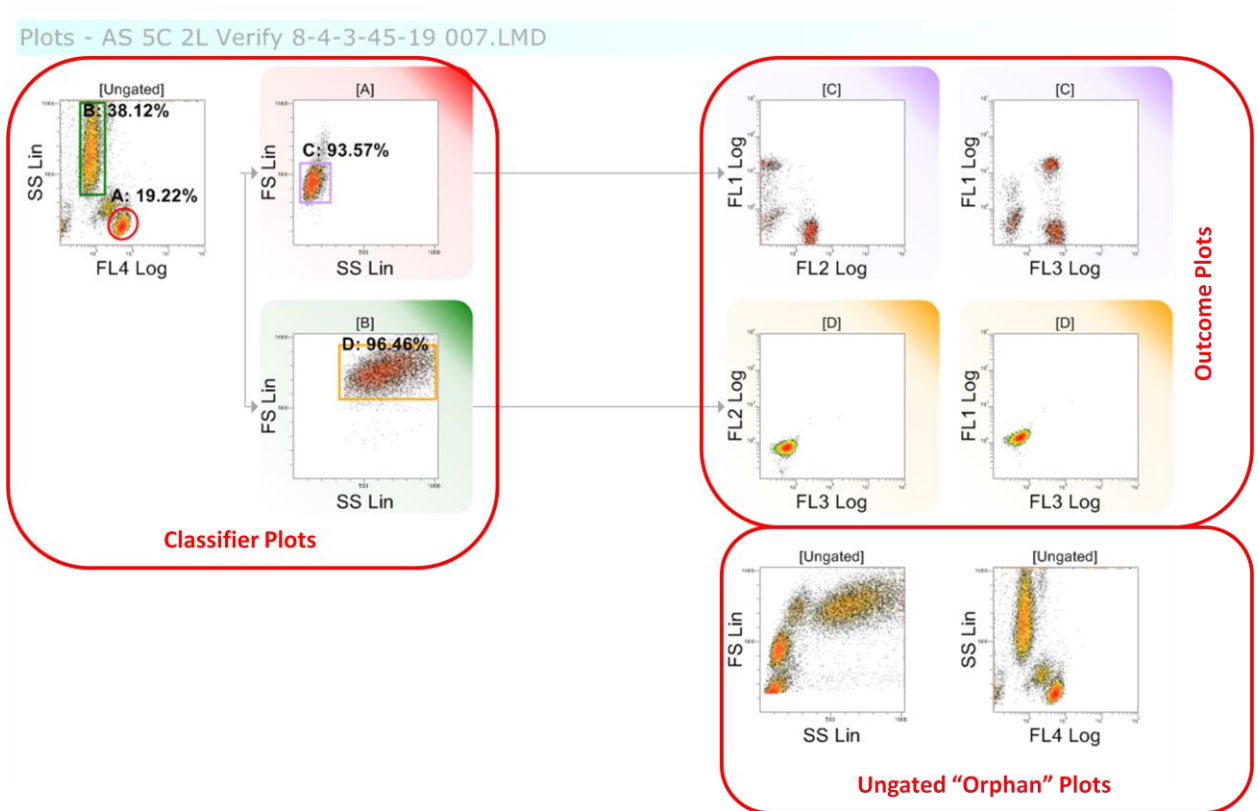

Example of plots on hierarchy format with only one classifier plot:

