



VenturiOne Ribbon

The **Ribbon** is the area at the top of the **Application Window** that presents an easy-to-browse main set of commands for VenturiOne software. The majority of commands are accessible from the **Ribbon**.

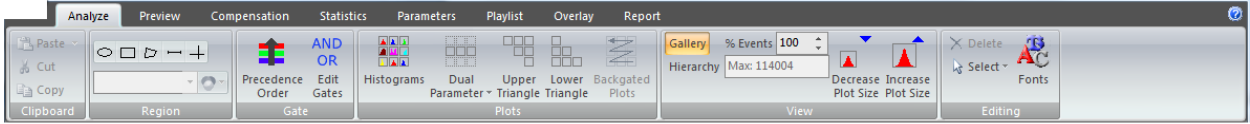


The **Ribbon** is organized into a set of tabs which represent the main functionality groups of the software:


Analyze tab – provides controls for analyzing and modifying the plots in the **Plots** area.

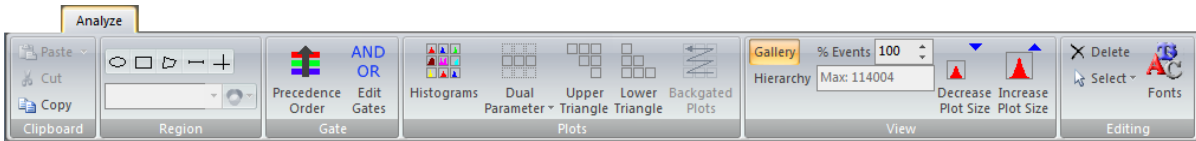


Preview tab - provides controls for the modification of the plots displayed in the **Previews** area.

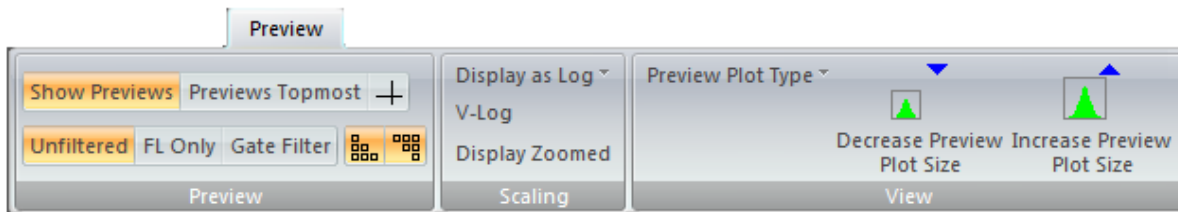

Compensation tab - provides controls allowing you to calculate fluorescence compensation.

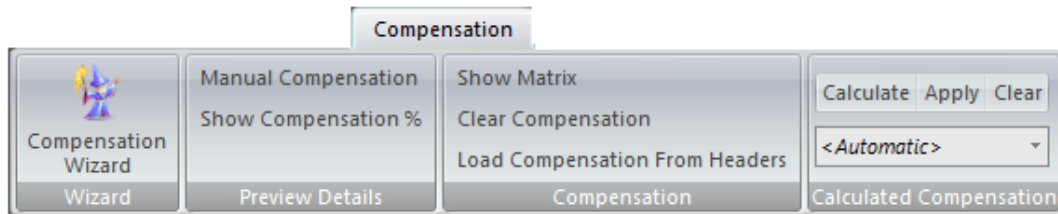


Statistics tab - provides controls for selecting which statistic is displayed on the plot(s) and which statistic(s) are presented in the associated results table.

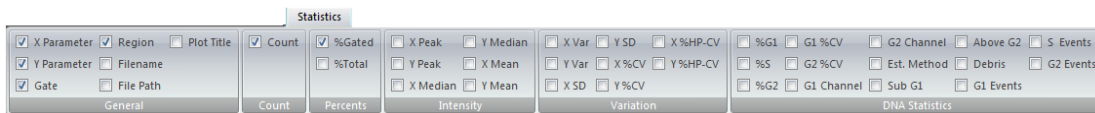


Parameters tab - provides controls for selecting the parameters displayed within the plots in the **Previews** area.




Playlist tab - provides controls for modification of all aspects of a **Playlist** including playlist items, LMD files, workspace files and compensation files.

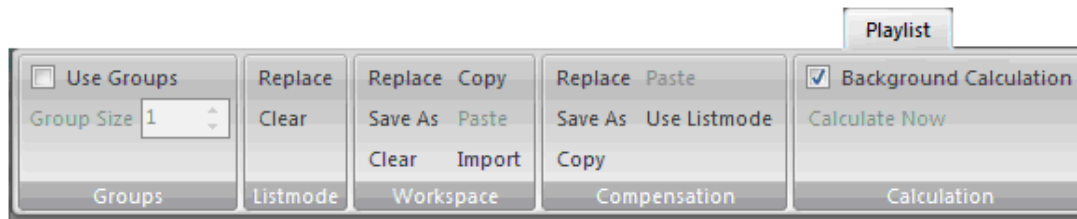


Overlay tab – provides controls for creating single and dual parameter plot overlay.

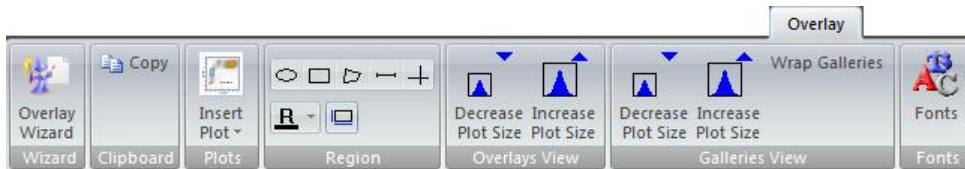

Report tab – provides controls for creating page based reports from a single playlist item or a group of playlist items.

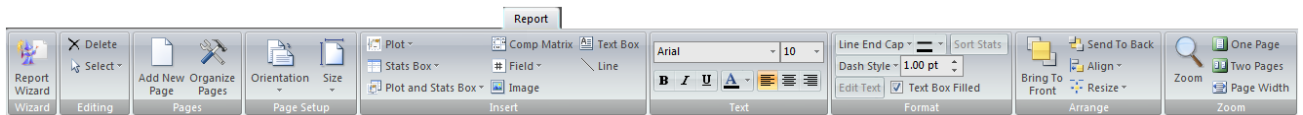

Plot Format tab - options here allow you to adjust settings on a per plot basis.

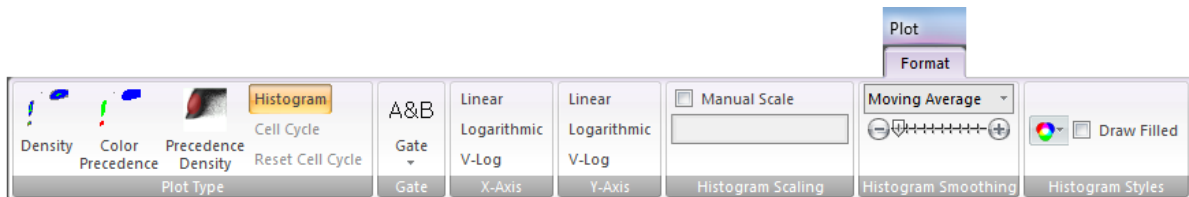

Region Format tab - The options here allow you to modify the appearance and type of region displayed on the selected plot. The region label can also be re-positioned and the regions can be tracked to other regions.

